


TONOVET Pro



TONOVET Pro tonometer provides **200° of positional freedom**


- **Downwards** measurement is now also possible
- **Probe release control** prevents accidental probe drops – even when the power is off

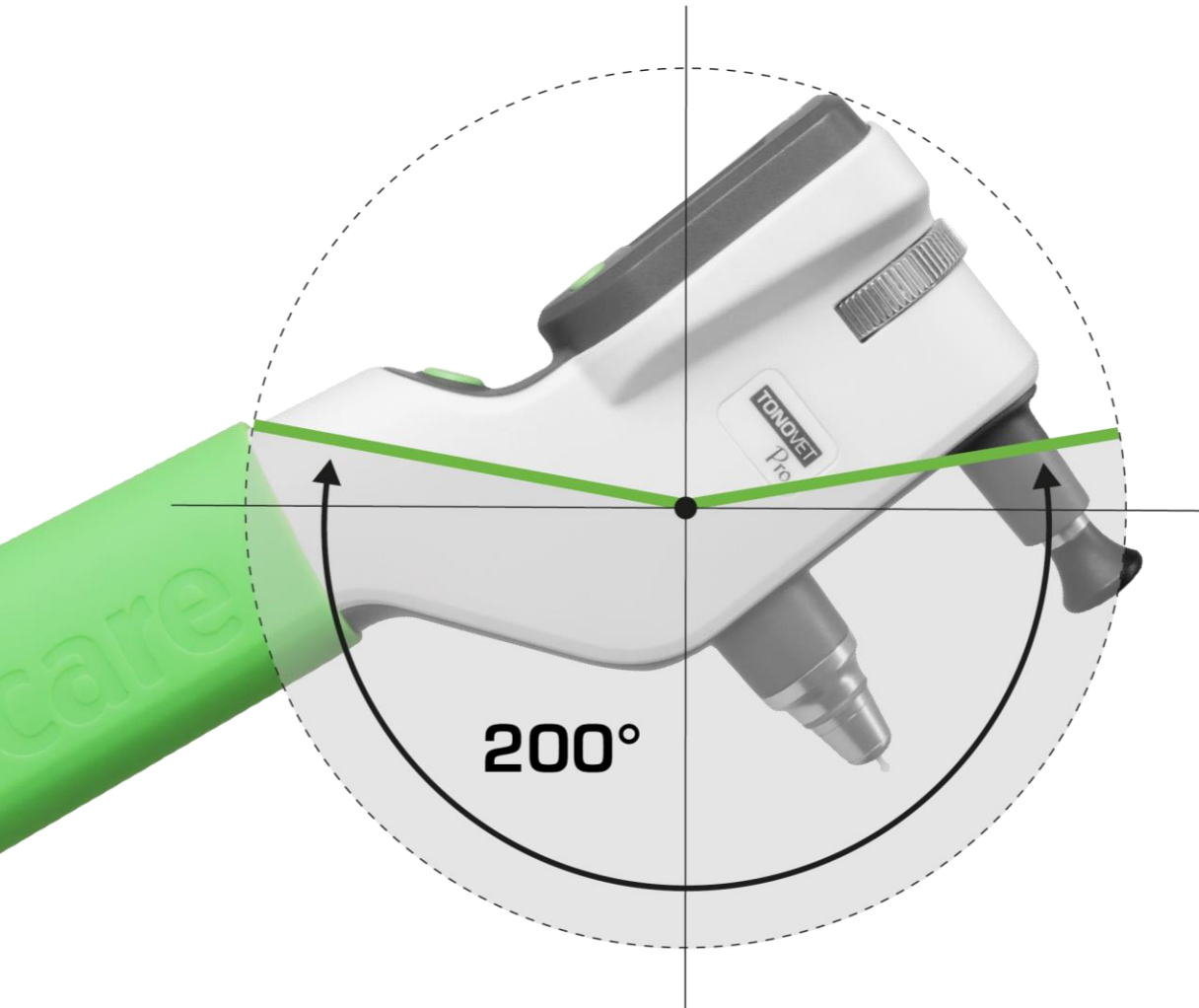

New Features

- 200° freedom of measurement
- Default and Quick Measure
- Improved probe reactivity
- Probe release control
- Intuitive user interface
- Alignment and position assistance
- Wireless printing via Bluetooth

TONOVET Pro

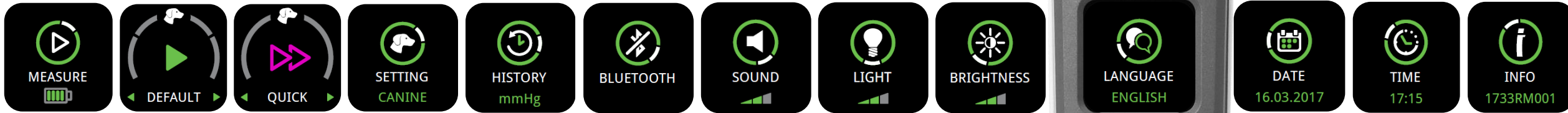

Intuitive User Interface

- Two navigation buttons + Select button
- Adjustable sound, probe base light and display brightness levels
- 12 languages


The **Measurement Button** is reserved for measurement only!

- You can access the measurement sequence anywhere by pressing the measurement button
- Series measurement mode



TONOVET Pro

- The TONOVET Pro tonometer provides faster IOP measurements
- Improved probe reactiveness
 - No delay in the probe launch
 - Probe speed is the same as in other products

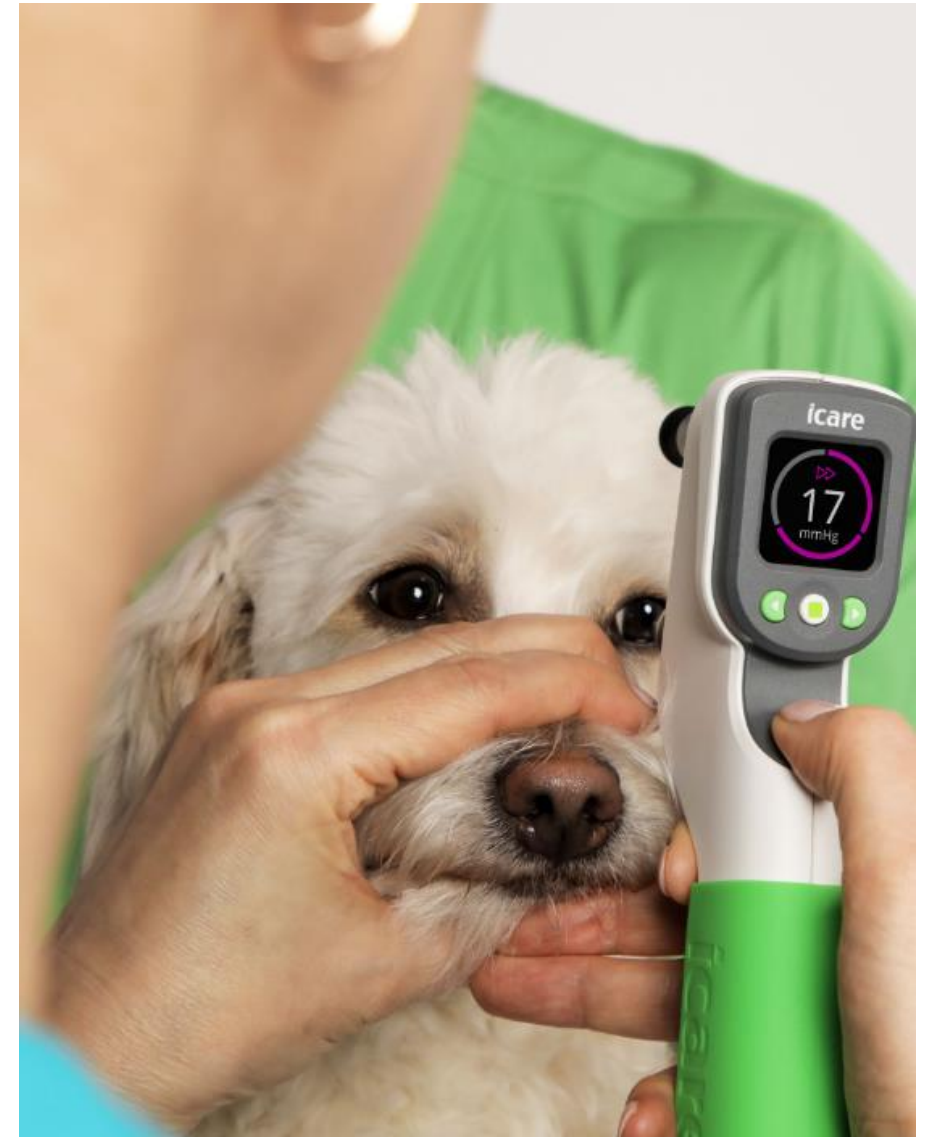


TONOVET Pro Features

The TONOVET Pro tonometer provides **two measurement options**

Default Measure

- 6 single rebound measurements
- Alignment assistant is ON
- **Highest accuracy**

Preferred for cooperative patients

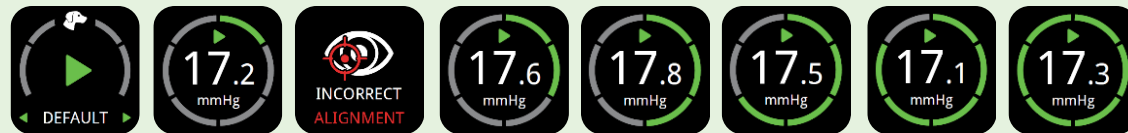

Quick Measure

- 2-3 single rebound measurements
- Alignment assistant is OFF
- **Easier and faster**

Suitable for challenging/uncooperative patients and inexperienced users

